

## Bench Mounted Laser Detection Head VIS



Bench Mounted Laser Detection Head VIS, #55-297

Stock **#55-297** **4 In Stock**

⊖ 1 ⊕ **A\$404<sup>00</sup>**

**ADD TO CART**

| Volume Pricing |                               |
|----------------|-------------------------------|
| Qty 1-5        | <b>A\$404.80</b> each         |
| Qty 6-24       | <b>A\$384.00</b> each         |
| Need More?     | <a href="#">Request Quote</a> |

### Product Downloads

### General

Detection Head **Type:**

**Typical Applications:**  
Ar-Ion, HeNe, HeCd, Nd:YAG, etc.

### Physical & Mechanical Properties

35.00 **Active Area Diameter (mm):**

70.00 **Outer Diameter (mm):**

8.00 Thickness (mm):

## Optical Properties

VS Wavelength:

Orange/Red (655nm), Broadband (600 - 730nm) Emission Color:

Band 1: 400 - 640nm  
Band 2: 800 - 1700nm Stimulation Range:

2 kW/cm<sup>2</sup> @ 1064nm, 7ns, 10Hz Minimum Stimulation, Pulsed:

## Electrical

Visible: 0.5 - 3 s (dependent on ambient light)  
IR: <0.5 s Persistence (Stimulation Removed):

<1 nW/cm<sup>2</sup> @ 450nm  
<25 μW/cm<sup>2</sup> @ 950nm Minimum Stimulation, Continuous:

100 W/cm<sup>2</sup> @ 512nm Maximum Stimulation, Continuous:

850 MW/cm<sup>2</sup> @ 337nm, 4ns  
60 MW/cm<sup>2</sup> @ 1064nm, 7ns Maximum Stimulation, Single Pulse:

## Regulatory Compliance

Compliant Reach 191:

Compliant RoHS 2015:

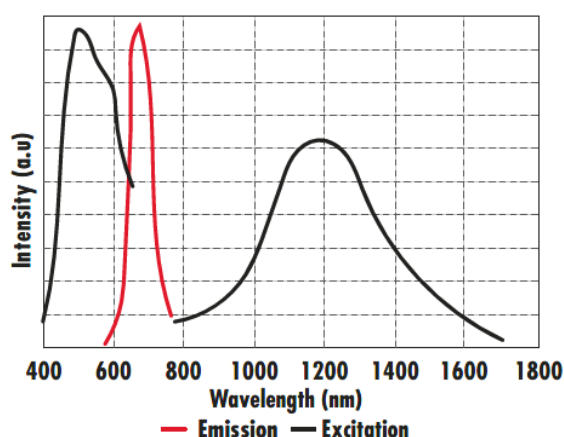
View Certificate of Conformance:

## Product Details

- Full Spectrum Coverage: UV, VS, IR Series
- 3 Mounted Formats Have Safe, Non-reflective Encapsulation
- Unique, No Pre-charge for IR Detection and No Fading During Use
- Flexibility for Either Transmission or Reflective Viewing

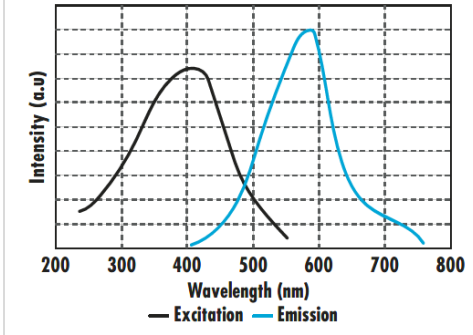
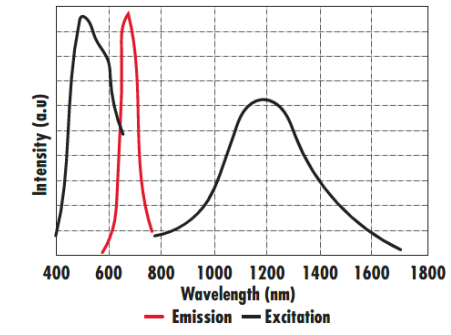
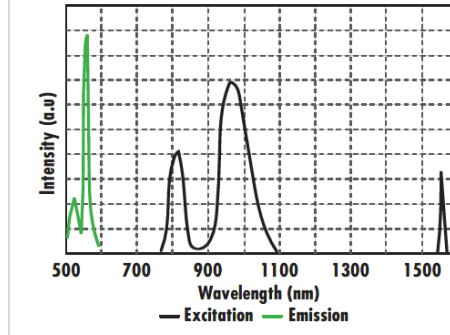
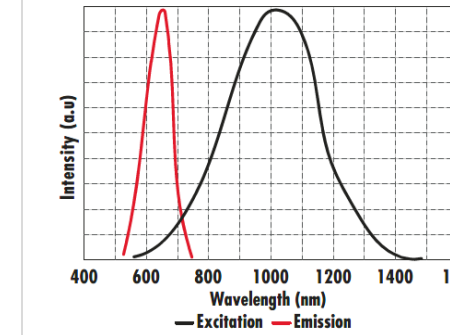
Laser Detection Products offer UV, visible, and IR laser users' greater performance and safety. They reduce problems associated with beam visualization, profiling, and alignment in many applications. Each range is available in three formats. Laser Detection Products' laminated credit card style is for low-power sources and reflective viewing only. The 25mm disk and clip-on wand style is used when frequent component positioning is required. The removable disk is positional at an optics location to enable precise alignment, while the wand format permits handling into the beam path. The optical bench-mountable head format has a large active area and ¼-20 threaded mounting for standard English post/post holder integration.

## Technical Information



VIS Detection Products

| Laser Detection Products |                                    |   |   |                     |
|--------------------------|------------------------------------|---|---|---------------------|
|                          | UV                                 | VS  | IR  | NIR                 |
| Stimulation Range        | 250 - 550nm                        | Band 1: 400 - 640nm<br>Band 2: 800 - 1700nm | Band 1: 790 - 840nm<br>Band 2: 870 - 1070nm<br>Band 3: 1550nm                           | 700 - 1400nm        |
| Typical Applications     | HeCd, Ar-Ion, tripled Nd:YAG, etc. | Ar-Ion, HeNe, HeCd, Nd:YAG, etc.            | 808nm, 820nm, 830nm, 880nm, 960 - 980nm Laser Diodes, Nd:YAG, 1550nm telecommunications | Nd:YAG, Fiber Laser |

|                                    |  |  |  |  |
|------------------------------------|--|--|--|--|
| Emission Color                     | Yellow (580nm), Broadband (490nm - 700nm)  | Orange/Red (655nm), Broadband (600 - 730nm)  | Green (550nm), other peaks at Red (673nm) and Blue (400nm)   | Orange/Red (655nm)   |
| Persistence (Stimulation Removed)  | 6 s - 4 mins (dependent on ambient light)  | <b>Visible:</b> 0.5 - 3 s (dependent on ambient light) <b>IR:</b> <0.5 s   | 800µs  | <50 ms   |
| Continuous (Minimum Stimulation)*  | <1nW/cm <sup>2</sup> @ 450nm & 365nm   | <1nW/cm <sup>2</sup> @ 450nm<br><25µW/cm <sup>2</sup> @ 950nm  | <2µW/cm <sup>2</sup> @ 808nm<br><175 nW/cm <sup>2</sup> @ 960nm<br><100µW/cm <sup>2</sup> @ 1550nm | 8µW/cm <sup>2</sup> @ 1064nm   |
| Pulsed (Minimum Stimulation)*      | <8W/cm <sup>2</sup> @ 337nm, 4ns, 20Hz <40W/cm <sup>2</sup> @ 337nm, 4ns, 1Hz  | 2 kW/cm <sup>2</sup> @ 1064nm, 7ns, 10Hz   | 250 kW/cm <sup>2</sup> @ 1064nm, 7ns, 10Hz   | N/A  |
| Continuous (Maximum Stimulation)   | 100W/cm <sup>2</sup> @ 512nm (all formats)   | 100W/cm <sup>2</sup> @ 512nm (all formats)   | 100W/cm <sup>2</sup> (all formats)   | 100W/cm <sup>2</sup> @ 1064nm (estimated)  |
| Single Pulse (Maximum Stimulation) | 130MW/cm <sup>2</sup> @ 337nm, 4ns (card only)<br>850MW/cm <sup>2</sup> @ 337nm, 4ns (other formats)<br>60MW/cm <sup>2</sup> @ 1064nm, 7ns (all formats) | 130MW/cm <sup>2</sup> @ 337nm, 4ns (card only)<br>850MW/cm <sup>2</sup> @ 337nm, 4ns (other formats)<br>60MW/cm <sup>2</sup> @ 1064nm, 7ns (all formats) | 35MW/cm <sup>2</sup> @ 1064nm, 7ns (all formats)   | 35MW/cm <sup>2</sup> @ 1064nm, 7ns (estimated)                                       |
|                                    |    |   |                |  |

\*Measured in darkened conditions