

[See all 18 Products in Family](#)

InfiniTube for Zeiss Objectives

See More by [Infinity Photo-Optical Company](#)



Standard (with no In-Line Attachment), InfiniTube, #56-125



Stock **#58-538** [CONTACT US](#)

-

1

+

A\$2,080⁰⁰

ADD TO CART

Volume Pricing	
Qty 1+	A\$2,080.00 each
Need More?	Request Quote

Product Downloads

General	
InfiniTube™ Z	Model Number:
Lens Accessory	Type:
Threading & Mounting	
C-Mount	Mount:

Regulatory Compliance	
-----------------------	--

RoHS 2015:

Reach 224:

Certificate of Conformance:

- **Mitutoyo** and **Achrovid™** Objectives Sold Separately

This in-line assembly unit is used to mount **Mitutoyo** or **Achrovid™** long working-distance objectives to C-mount cameras, yielding magnifications not possible with ordinary video lenses, such as 50X magnification. The InfiniTube™ system can be configured in a variety of fashions. It can be utilized alone (InfiniTube™ Standard **#56-125**) for applications where illumination is provided from another source such as a backlight or directional illumination.

The first in-line attachment ([#56-126](#)) has been designed for high performance with objectives having a magnification of 10X or below and is not suggested for use with higher magnification objectives. The second in-line unit ([#56-191](#)) has been designed to work with increased efficiency when used with magnifications of 10X or higher, but can be utilized with all magnifications. An InfiniTube™ In-line Assembly ([#54-590](#)) is also offered, comprised of the InfiniTube™ Standard and the higher magnification in-line attachment. This Assembly can be used with all magnifications. A filter holder, which attaches to the InfiniTube™ itself, can be added to the systems to introduce 25mm diameter filters. These lenses are commonly used in microscopy and wafer inspection.


By utilizing an in-path beamsplitter and side-port illumination tube with collimating optics, the In-Line Assembly (#54-590) provides incident illumination parallel to the optical axis (coaxial). Coaxial illumination requires an optional [Fiber Optic Illuminator](#) and a 1/4" diameter [Fiber Optic Light Guide](#).

A version of the InfiniTube™ is specifically designed only for Zeiss ICS-corrected objectives (#58-538, 160mm tube lens design). Shorter than the Standard model, this model has color-compatible correction for use exclusively with Zeiss ICS objectives. The RMS adapter (#53-786) is required to connect #58-538 to the Zeiss EC Plan-Neofluar® objectives.

The RMS adapter ([#53-786](#)) is required to connect [#58-538](#) to the Zeiss ICS objectives. Since even the slightest vibration will affect such small fields-of-view seen with the InfiniTube™ we recommend [vibration isolation equipment](#) when using objectives greater than 10X in magnification. We also recommend the use of a focusing mechanism (such as the Rack and Pinion Movement mount [#54-794](#)) when using the lens alone.

Mounting clamp (#57-788) is included with #54-590, #56-125 and #58-538. Note that an objective, objective adapter, light guide, and illuminator are not included.

Technical Information

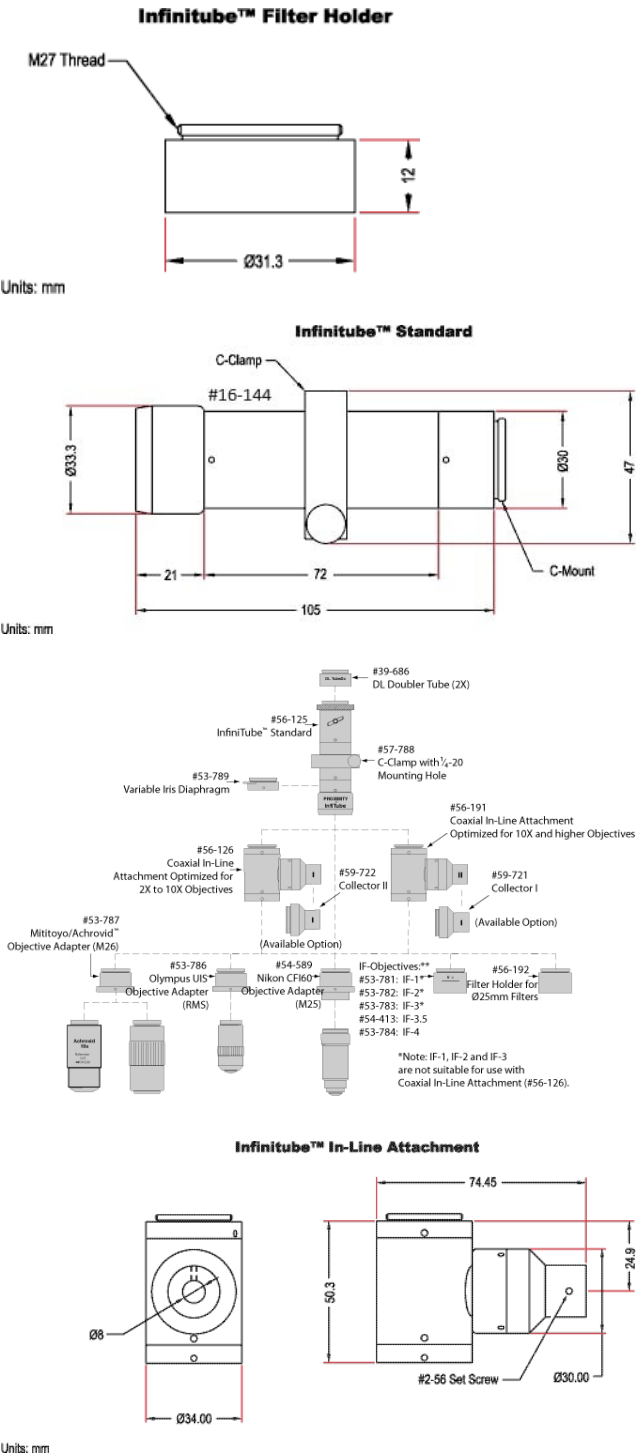

#54-562	
Item	Stock No.
InfiniTube In-Line Assembly Unit	#54-590

Hitachi KP-D20B Camera	#55-837
15.1" LCD Monitor (US Version Only)	#56-522
Y-C Cable (6 ft.)	#39-255
Mitutoyo/Achrovid™ Objective Adapter	#53-787
MI-150 Fiber Optic Illuminator	#55-718
Fiber Optic Light Guide (48" L)	#39-367
Fiber Optic Adapter	#38-944
Standard Boom Stand	#54-120
Rack & Pinion Mounting Plate	#54-123