


[See all 17 Products in Family](#)

InGaAs FC Fiber Balanced Photoreceiver, 900-1700nm

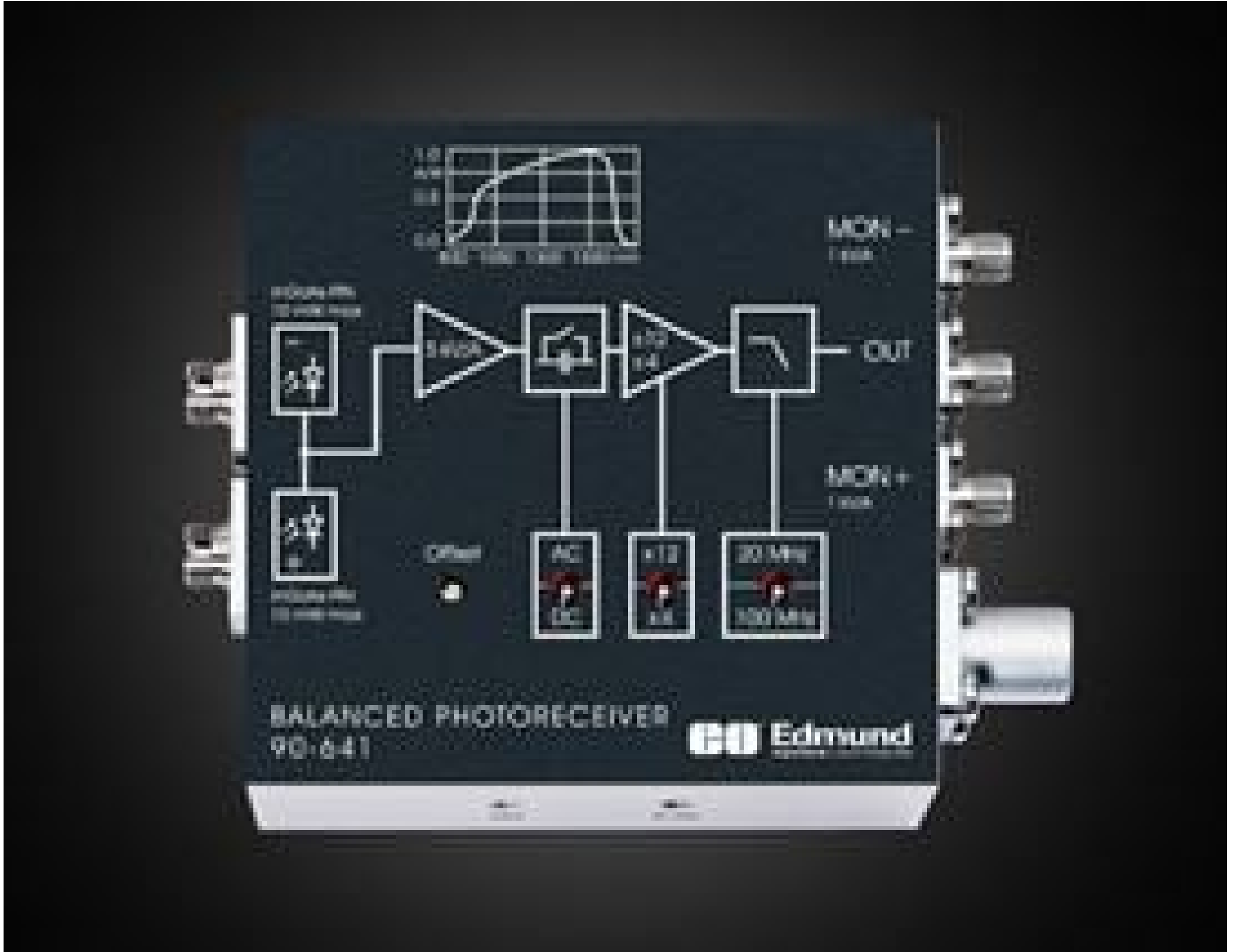

Stock #90-641 **NEW** 1 In Stock

1 **A\$5,704⁰⁰**

ADD TO CART

Volume Pricing

Qty 1+	A\$5,704.00 each
Need More?	Request Quote

Product Downloads

General

Note:

Includes:
Lemo®3-pin connector
(3) Adapter SMA (male) to BNC (female)
Datasheet

Physical & Mechanical Properties

Weight (g):

350

Dimensions (mm):

Case Size: 80 x 80 x 30.5

Optical Properties

Spectral Range:
900 - 1700 nm

Sensor

Detector Type:
InGaAs PIN photodiode

Electrical

Transimpedance Gain (Ω):
 2×10^4 or 6×10^4 (switchable)

Noise Equivalent Power NEP ($W/Hz^{1/2}$):
 4×10^{-12} @880nm

Bandwidth (-3 db):
100 MHz

Output Signal:
 ± 1.0 V at 50 Ω load (for linear gain and low harmonic distortion), maximum ± 2.0 V at 50 Ω load

Conversion Gain (V/W):
 20×10^3 or 57×10^3 (switchable)

Common Mode Rejection (dB):
55

Hardware & Interface Connectivity

Power Requirement:
 ± 15 V (± 14.5 V... ± 16.5 V) -90 / +120 mA

Power Supply:
Power Supply Required and Sold Separately.
USA: [#59-180](#)
Europe: [#59-180](#)
Japan: Not Available
Korea: Not Available
China: [#59-180](#)

Environmental & Durability Factors

Operating Temperature ($^{\circ}C$):
0 to +60

Regulatory Compliance

RoHS 2015:
[Compliant](#)

Certificate of Conformance:
[View](#)

Product Details

- Subtracts Two Photodiode Signals for Differential Detection
- Improved Signal to Noise Ratio (SNR) for Weak or Modulated Optical Signals
- High Common-Mode Noise Suppression for Improved Measurement Sensitivity and Accuracy
- Available in Si and InGaAs models for UV-VIS and NIR spectral ranges

Balanced Photoreceivers use true differential detection by subtracting the photocurrents from two matched photodiodes, producing a single electrical output proportional to the difference in optical power between the two inputs. This suppresses common-mode noise, such as laser intensity fluctuations, improving SNR and measurement sensitivity. Balanced Photoreceivers are engineered with a low-noise transimpedance amplifier, which ensures stable, consistent performance for precision optical measurements. Available in Si and InGaAs models for UV-VIS (320-1000nm) and NIR (800-1700nm) spectral ranges, these photoreceivers are ideal for coherent optical detection, interferometry, spectroscopy, and optical coherence tomography (OCT).

Note: Power supply sold separately. Please see specifications for more details.